

RADD5 - Detoxification, Relaxation and Rejuvenation Tool






Operation Guide


Disclaimer RADD5 is a consumer relaxation product, not a substitute for professional medical care. It does not intend to diagnose, treat cure or prevent any medical conditions or diseases.

Warning DON'T use RADD5 if you have pacemaker or artificial heart.

During the initial use of the RADD5, each individual may experience different reactions such as fatigue, dizziness or muscle ache, etc. These are the signs of detoxification. Normally, the reactions will last for 1 to 2 days (some might take longer).

Please follow the steps below to operate the RADD5:
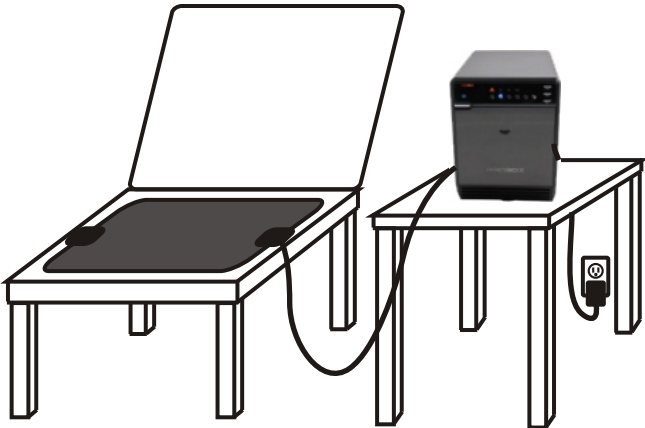
<p>A. </p> <p>B. </p> <p>C. </p> <p>D. </p>	<p>Step 1. Make sure that you have the following items inside the packaging box:</p> <ul style="list-style-type: none">A. RADD5 UnitB. Charcoal Colored MatC. Ion CableD. Power Adaptor
	<p>Step 2. Place the charcoal colored mat on a chair or bed, and the RADD5 unit on a table next to the chair or bed. And locate an electrical power outlet which is close to the setting.</p> <p>Step 3. You may sit or lay on the mat with your cloth on when using the RADD5. To get better result, wear cotton cloth and drink at lease one cup of water before and one cup of water after using. After turn on the Ion Output, close your eyes and feel the subtle energy flow inside your body. Sometimes, it may take few minutes to feel that.</p>

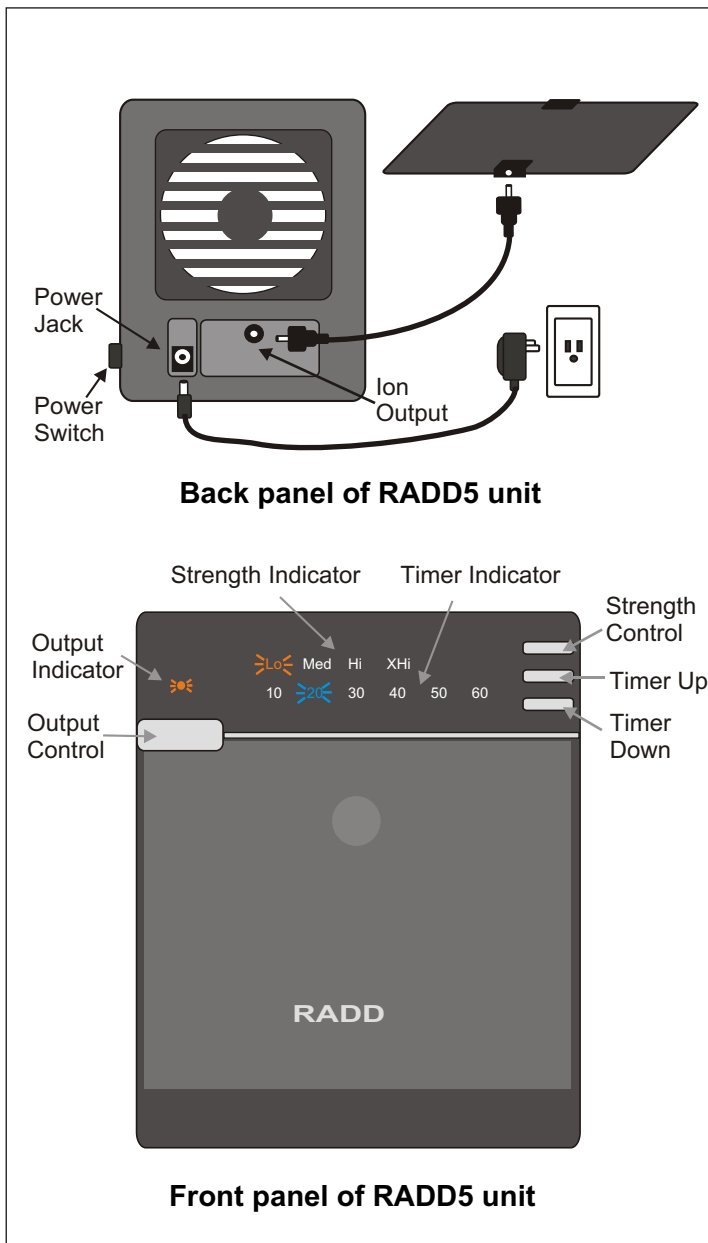

Step 4. Plug one end of the Ion Cable to the Ion Output and the other end to the charcoal colored mat.

Step 5. Plug the Power Adaptor to an electrical outlet and the other end to the Power Jack. Then turn the Power Switch on. The Timer Indicator 20 and the Strength Indicator Lo will light on. And you will hear four beeps. The RADD5 now is in Default setting.

Step 6. To use Default setup (Ion Output Strength Lo and Ion Output on for 20 minutes), go to step 9.

Step 7. Set the Ion Output Strength to higher level by pressing the Strength Control button: once for Medium, twice for High, three times for Extra High. Press again back to Low.

Step 8. Set the Timer by pressing the Timer Up button or the Timer Down button to the desired number of minutes.

Step 9. Press the Output Control button once to turn on the Ion Output. The Output Indicator will illuminate red when the Ion Output is on.

Step 10. Press the Output Control button again to turn off the Ion Output, or it will turn off automatically when the timer is up.

Step 11. Turn the power Switch off when not using the RADD 5.

Please start to use the RADD5 at LOW Output once or twice a day. After a week or two, if you don't feel any uncomfortable, switch to Higher Output level. To use it before bedtime will improve your sleeping quality, after exercise will relieve the tight muscle, when under stress will reduce your body tension, improve mood and ease anxiety. If you feel strongly uncomfortable when using it, stop using and contact your vendor.